



Experience the
joy of walking


Fillauer
COMPANIES

K2 FOOT

Fillauer proudly presents its newest foot uniquely designed for the functional level 2 community ambulator. The Fillauer K2 Foot provides a flexible keel and multi-axis design assembled in a removable foot shell. Let your patients experience the joy of walking!



K2 FOOT

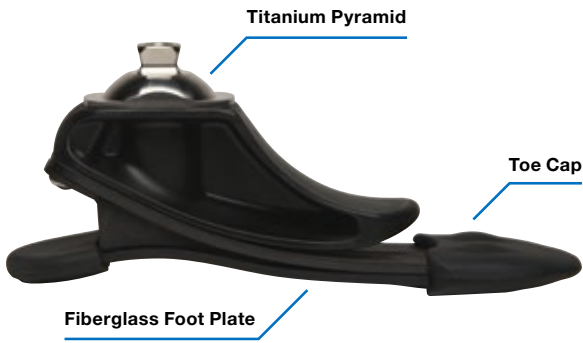

The new K2 Foot from Fillauer is engineered from the ground up for the functional level 2 patient. While other K2 designs typically fall into two categories, stiffer SACH feet and softer K3 designs, the Fillauer K2 Foot is a new design specifically for the functional level 2 patient. Patient testing shows K2 level patients walk with a more controlled gait pattern and need more stability at midstance. With this in mind, we engineered the K2 Foot with a unique flexible keel to provide better control and stability. When the patient is ready to go, so is the K2 Foot with a compliant forefoot that loads easily and provides support through toe off.

PRODUCT FEATURES

- Specifically engineered for K2 patients
- Unique flexible keel/multi-axis design
- Flexible, impact reducing design
- Roll-Over Shape (ROS) designed for a natural gait
- Compliant dynamics for comfort and stability
- Positive lock foot shell interface for stability
- Simple, maintenance-free design
- Foot comes assembled in foot shell for an easy, “out-of-the-box” fit

SPECIFICATIONS

- Patient weight: Up to 275 lbs. (125 kg)
- Foot weight (26 cm): 14 oz. (395 g)
- Build height: 3.5 in. (9 cm)
- 2 year warranty/30 day patient trial
- Functional level: K2

L-CODES*

- L5972** All lower extremity prostheses, flexible keel foot
- L5986** All lower extremity prostheses, multi-axial rotation unit (“MCP” or equal)

***Suggested L-Codes are provided as a reference only. It is the responsibility of the practitioner and/or patient care facility submitting this claim to determine and submit the correct coding within all established rules and guidelines.**

WARRANTY

- 24 Months from date of patient fitting
- Foot Shell—6 months from date of patient fitting

SATISFACTION GUARANTEE

- 30 days from date of patient fitting

ORDERING INFORMATION

134-CCC-XX-X K2 Foot
Example: 134-13L-25-B = Color 13, Left, 25 cm, 160 lbs.

To order, choose the color (CC 03, 09, or 13). Then, select the side (left or right).

	CC 03	CC 09	CC 13
Left	03L	09L	13L
Right	03R	09R	13R



		22 cm	23 cm	24 cm	25 cm	26 cm	27 cm	28 cm	29 cm	30 cm
99-130 lbs.	45-59 kg	22-A	23-A	24-A	25-B	26-B	27-B	28-B	29-D	30-D
131-164 lbs.	60-74 kg	22-A	23-A	24-A	25-B	26-B	27-C	28-C	29-D	30-D
165-199 lbs.	75-90 kg	22-B	23-B	24-B	25-C	26-C	27-C	28-C	29-D	30-D
200-235 lbs.	91-107 kg	22-B	23-B	24-B	25-C	26-C	27-D	28-D	29-E	30-E
236-275 lbs.	108-125 kg	22-C	23-C	24-C	25-D	26-D	27-D	28-D	29-E	30-E

K2 Foot includes foot shell. The foot model recommendations on the selection chart are based on average patients engaged in normal activities. If your patient is more active than average or lifts heavy loads as part of a daily routine, then selecting the next weight category higher is advised. If you are not sure which foot to select, please consult with one of our technical service representatives.