

# SPL2 Clinical Trial Orthosis


from Fillauer LLC

Fillauer is pleased to introduce The SPL2 Clinical Trial Orthosis. This innovative product allows the user the opportunity to wear and see the Benefits of the Fillauer SPL2 Orthotic Knee Joint prior to purchase. The CTO is easily adjustable for a Wide Variety of body types and knee centers.

## FEATURES

- 4 Modes of Operation
- Easily Adjustable Knee Center Height parameters
- Available in Left & Right
- Ability for the Clinician and Patient Candidate to **'Try Before You Purchase'** to achieve Best possible Outcome

## FITTING KNEE CENTER PARAMETERS

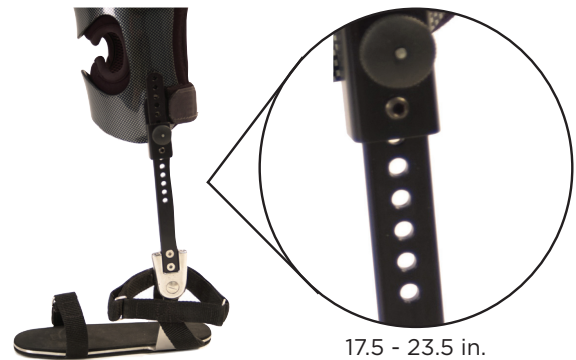

17.5 - 23.5 in.

## LOCKING MODES


All SPL2 Knee Joints are being executed with the new modified Satellite. The Satellite has the following features:

- Automatic Locking (Position 1)
- Manual Unlock (Position 2)
- Free Motion Mode (New) (Position 3)
- Manual Lock (Position 4)

NOTE: If the mode switch is pressed when in Free Motion Mode, the satellite will drop into the Lock Mode. The Automatic Lock must be re-engaged with the slide lever.

## ORDERING INFORMATION

- 023902** SPL2 Clinical Trial Orthosis, Left  
**023903** SPL2 Clinical Trial Orthosis, Right

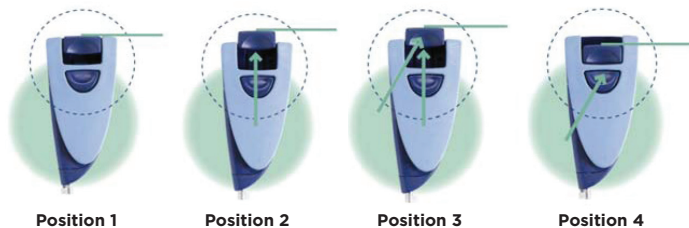

Position 1

Position 2

Position 3

Position 4

