



TASKA[®]



The world's only rugged,
water-resistant, multi-
articulating hand.



FLEXIBLE FINGERS

The fingers on the Taska hand are flexible and can laterally spread to allow the hand to effectively grip a wide variety of objects firmly.

CONFIDENT GRIP

Hold things with more confidence. Soft pads between the fingers and on fingertips, complement the flexible fingers to provide unparalleled grip security.

MOTORIZED THUMB ROTATION

High speed thumb rotation mimics the motion of the natural hand and allows nearly instant change of grip pattern.

THE TASKA® LOW PROFILE WRIST OPTION

A low profile wrist option extends water resistance up the forearm. The LP Wrist also allows for 90° of passive rotation.

KNUCKLE BREAKAWAY


Knuckles break away under large loads to prevent damage and are user resettable to provide additional robustness.

GRIP CYCLE BUTTONS

Simple touch-and-go, easy-to-use grip cycle buttons provide access to grips. EMG grip patterns can also be disabled at the touch of a button.

INTEGRATED FLEXIBLE WRIST

Every model of the Taska hand has a built-in, flexible wrist with three easily lockable positions, as well as center-biased free flexion.

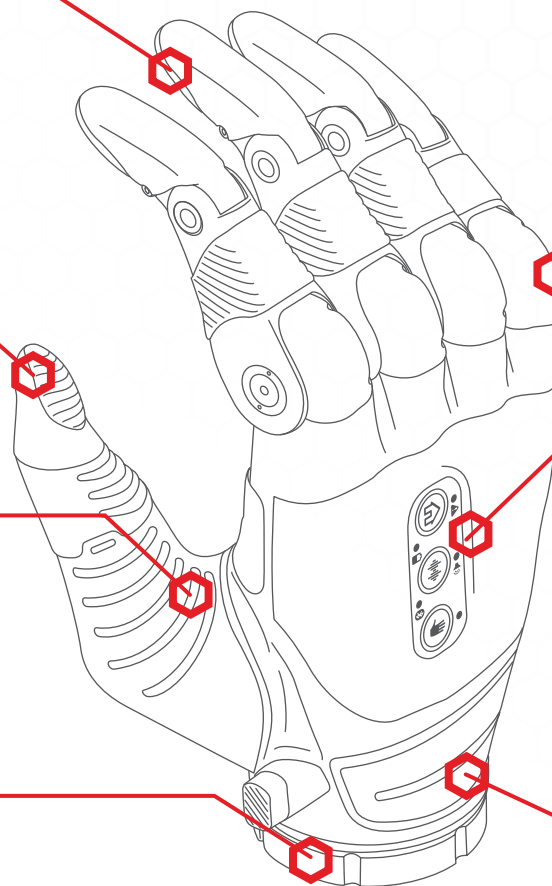

8 1/4



Black



Sand



White



7 3/4

Fillauer
Motion Control

801.326.3434

Exclusive Distributor for U.S and Canada

www.fillauer.com/taska